


DEFIGARD 400

Portable Biphasic Defibrillator

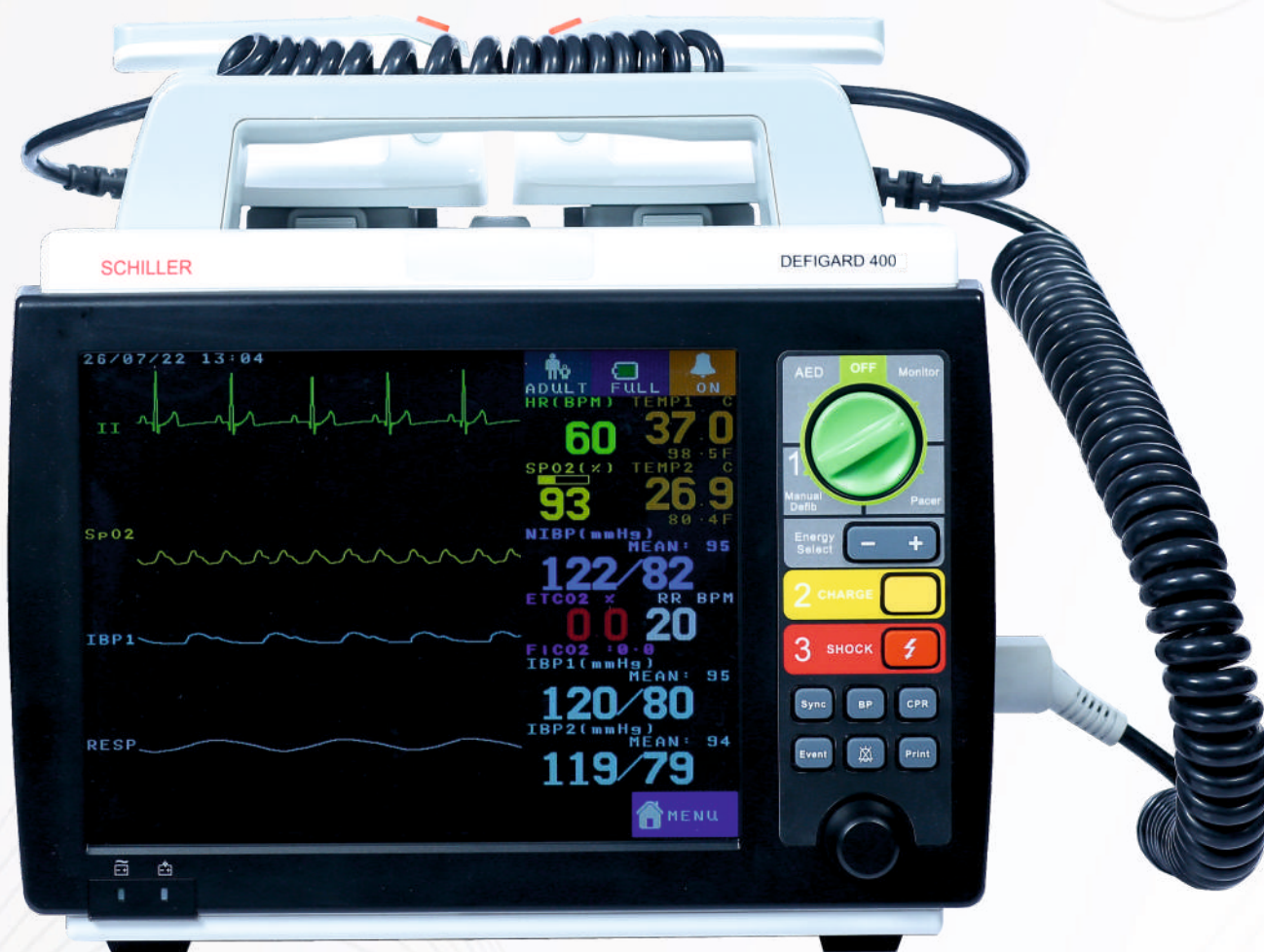

SCHILLER

The Art of Diagnostics

DEFIGARD 400

DEFIGARD 400 is a versatile biphasic defibrillator. Specifically designed for critical care, it combines portability, precision and power.



- ❖ The DEFIGARD 400 is an advanced portable defibrillator with an integrated monitor, designed to meet the needs of healthcare professionals in both in-hospital and pre-hospital environments.
- ❖ Built with reliability, ease of use, and innovative features, it supports rapid and effective responses in critical situations.



ENHANCED USABILITY AND PERFORMANCE



Wide Display: 10" color TFT display with LED backlight ensures clear visibility in any condition




User-Friendly Interface: Direct function keys and a rotary knob allow for easy mode selection and operation



Visual Instructions: On-screen guidelines provide step-by-step instructions, enhancing user confidence and accuracy



Memory Capacity: Internal SD Card: Stores up to 50 defib events and 50 nos. 12-lead ECGs. External USB: 1000 defib events

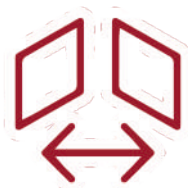

Extended Monitoring: 72 hours tabular and graphical trend monitoring



Performance: Unit operating time on fully charged battery: backup time 3 hours and no. of shocks more than 100

Energy Selection Through Paddles: The DEFIGARD 400 allows healthcare providers to select energy levels directly through the paddles, streamlining the workflow and ensuring swift responses during critical moments.


Contact Indicator: The paddles are equipped with two LEDs to show contact status, while the monitor also displays a "Contact OK" message for added assurance. This feature eliminates guesswork and ensures optimal energy delivery.

Child Mode Auto-Detection: Automatically switches to child mode when adult paddles are removed, simplifying paediatric care without the need for manual intervention.

EXTENDED FUNCTIONALITIES FOR VERSATILE USE

Enhance the DEFIGARD 400 with optional upgrades to tailor the device to specific clinical needs:

- ✦ **Arrhythmias:** Unit can be upgraded for arrhythmia detection from ECG.
- ✦ **AED Mode:** Simplifies operation during emergencies with automated external defibrillation.
- ✦ **External Pacemaker:** Offers advanced cardiac care with optional pacing functionality.
- ✦ **Internal Defibrillation:** Equipped with internal paddles for surgical settings.
- ✦ **Comprehensive Monitoring:** 360J, Includes 12-channel ECG, SpO₂ (Digital/ Nellcor), NIBP, Temp, Dual IBP, EtCO₂.
- ✦ **Advance CPR Feedback:** Ensures effective chest-compressions with Metronome - real time feedback on depth & Rate for CPR.
- ✦ **Integrated Printing:** Onboard thermal printer for immediate data documentation.

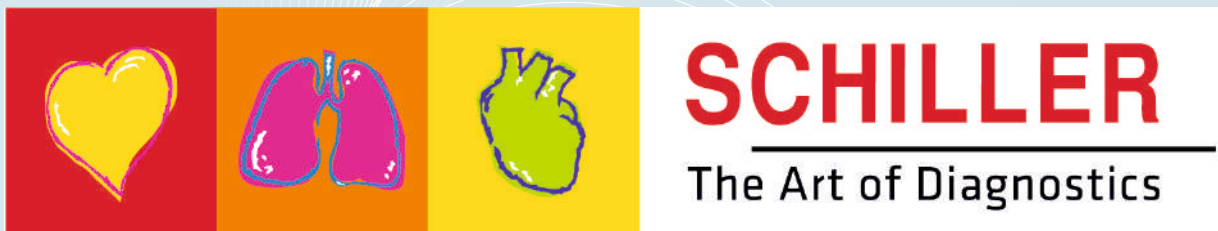

SCHILLER, founded in 1974 by Alfred E. Schiller is a successful group with 30 subsidiaries and a global sales network. Today, SCHILLER is a world-leading manufacturer and supplier of devices for cardiopulmonary diagnostics, defibrillation and patient monitoring as well as software solutions for the medical industry.

Over the past 25+ years, Schiller India has established itself as a leader in the medical technology sector. With 450+ employees, a state-of-the-art ISO 13485-certified production centre in Puducherry, and multiple R&D centres across the globe, Schiller India makes advanced healthcare equipment accessible through a network of 100 sales and service dealers across more than 45 locations. Our product range includes Critical Care, Anaesthesia, Emergency Care, Cardiology, Respiratory Diagnostics, Radiology, and Robotics.

Scan here
to know more:



Part no.: B1.100009  09/2025 Ver. no.: 1.3

All registered trademarks acknowledged.

SCHILLER
The Art of Diagnostics

Swiss H.Q.: SCHILLER AG, Altgasse 68, P. O. Box 1052, CH -6341 Baar, Switzerland

Indian H.O.: Schiller Healthcare India Pvt. Ltd., Advance House, Makwana Road, Andheri (East), Mumbai-59. Tel : + 91 - 22 - 61523333

Factory: No. 15/5 & 15/6, Vazhuthavur Road, Kurumbapet, Puducherry - 605 009, India

Regional Offices: Mumbai: 022-61523333 Email: sales.west@schillerindia.com Delhi: 011-41062067 / 09312432205 Email: sales.north@schillerindia.com Kolkata: 033-23593102 / 033-23593103 Email: sales.east@schillerindia.com Chennai: 044-28232648 / 044-28311021 Email: sales.south@schillerindia.com Ahmedabad: 079- 35337949 / 079- 35337950 Email: sales.west@schillerindia.com Puducherry: 0413- 220 2680 / 0413- 229 2940 Email: sales.south@schillerindia.com

CIN : U33110MH1997PTC111307

 sales@schillerindia.com, help@schillerindia.com  www.schillerindia.com  1800 209 8998

Follow us on :  /SchillerHealthcareIndia  /schiller_healthcare_india  /schillerhealthcareindia  /schiller_healthcare_india  @schillerindia1074