

CARDIOVIT CS-200 Office

Wireless ECG acquisition




CARDIOVIT CS-200 OFFICE

The CARDIOVIT CS-200 Office offers cutting-edge diagnostic capability for 12-lead resting and stress wireless ECG acquisition. The choice is yours: Whether only the software for any laptop or as a package integrated in a medical PC, it's the perfect solution for GPs, clinics and hospitals.



MS-12 BLUE

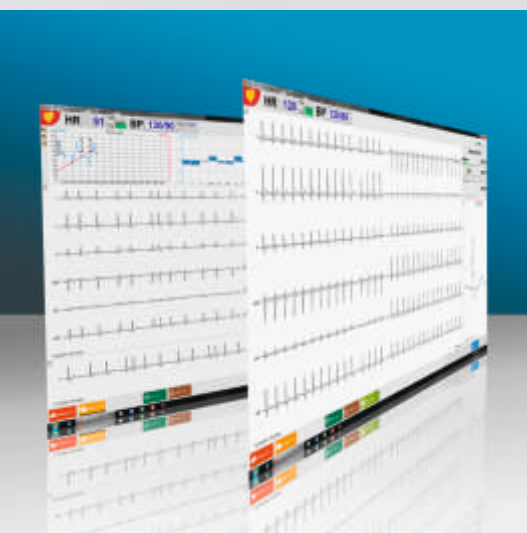

**SWITCH BETWEEN TWO CUSTOMISABLE
SCREEN LAYOUTS AT THE TOUCH OF A
BUTTON**

ECG ACQUISITION WITH MAXIMUM FREEDOM

With the introduction of the MS-12 blue, SCHILLER provides a new solution that greatly simplifies ECG acquisition.

- ✦ Allows to freely operate in any environment without worrying about cable length.
- ✦ Weighing only 115g and fitting in one hand, the MS-12 blue sits on the bed or on the patient body during ECG acquisition.
- ✦ Motion artefact is minimized without an analogue cable connection and an extremely lightweight patient cable.
- ✦ One battery charge guarantees continuous operation for more than 36 hours and more than a week with the power-saving mode.

ADVANCED CLINICAL TOOLS

The CS-200 Office can be supplemented with high-end software features:

- ✦ ETM interpretation
- ✦ ETM sport
- ✦ 12-lead exercise ECG
- ✦ 3D vector ECG
- ✦ Arrhythmia detection and classification
- ✦ Risk factors

USER-DEFINED WORKFLOW (WITH PRINTOUTS OR PAPERLESS)

The CARDIOVIT CS-200 Office fulfills the needs of every user thanks to a customizable user interface, stage and report printout with laser printer, and seamless connectivity to EMR, PACS, HIS or SCHILLER SEMA data management system.

SOLUTIONS FOR ANY ENVIRONMENT

The CS-200 Office can be easily installed on your laptop or PC, as well as purchased as a ready-to-use system including a robust trolley and a 21.5fi touch panel PC with pre-installed software (see picture).

Note:

CARDIOVIT CS-200 Office is also available with MS-12 USB (USB connection).

High-resolution and easy-to-clean
21.5fi touch screen.

Intuitive and clear user interface
suitable for touch screen.

Soft handle positioned at an
adjustable height.

Support for PC power
supply and cables.


Transparent bin to store and
transport medical items.

User-defined screen layouts
and workflow.

MS-12 blue: Small but
powerful wireless ECG
acquisition module with
on-screen signal view and
hookup advisor.

Easy-glide wheels to smoothly
maneuver the trolley.

*PC and Trolley are for representation purpose only.

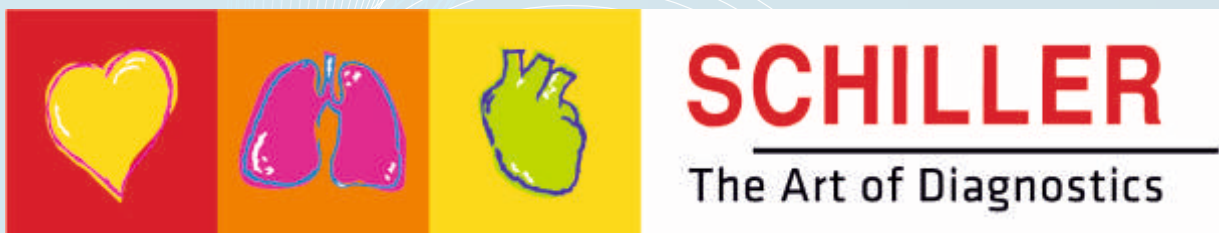

About SCHILLER

SCHILLER A. G., a company based in Switzerland, was founded by Alfred Schiller in 1974. Over the last 50 years, it has evolved into a successful Global Group with 31 subsidiaries and dealers across 30 countries. Today, SCHILLER is a world-leading manufacturer and supplier of devices for cardiopulmonary diagnostics, defibrillation and patient monitoring as well as software solutions for the medical industry.

Schiller in India was founded as a joint venture in 1997, and since then has established itself as a leader in the medical technology sector. With 450+ employees, a state-of-the-art ISO 13485-certified production centre in Puducherry, multiple R&D centres, and 17 offices across India; Schiller makes advanced healthcare equipment accessible through a network of 100 sales and service dealers across more than 45 locations. Our product range includes Critical Care, Anaesthesia, Emergency Care, Cardiology, Respiratory Diagnostics, Radiology, and Robotics.

Scan here
to know more:



Part no.: B3.100019 11/24 Ver. no.: 1

All registered trademarks acknowledged.

SCHILLER
The Art of Diagnostics

Swiss H.Q.: SCHILLER AG, Altgasse 68, P. O. Box 1052, CH -6341 Baar, Switzerland

Indian H.O.: Schiller Healthcare India Pvt. Ltd., Advance House, Makwana Road, Andheri (East), Mumbai-59. Tel.: + 91 - 22 - 61523333

Factory: No. 15/5 & 15/6, Vazhuthavur Road, Kurumbapet, Puducherry - 605 009, India

Regional Offices: Mumbai: 022-61523333 Email: sales.west@schillerindia.com Delhi: 011-41062067 / 09312432205 Email: sales.north@schillerindia.com Kolkata: 033-23593102 / 033-23593103 Email: sales.east@schillerindia.com Chennai: 044-28232648 / 044-28311021 Email: sales.south@schillerindia.com Ahmedabad: 079- 35337949 / 079- 35337950 Email: sales.west@schillerindia.com Puducherry: 0413- 220 2680 / 0413- 229 2940 Email: sales.south@schillerindia.com

CIN : U33110MH1997PTC111307

sales@schillerindia.com, help@schillerindia.com www.schillerindia.com 1800 209 8998

Follow us on : [f /SchillerHealthcareIndia](#) [i /schiller_healthcare_india](#) [in /schillerhealthcareindia](#) [@schiller_healthcare_india](#) [@schillerindia1074](#)