


CARDIOVIT AT-102 G2

Easy, smart and state-of-the-art

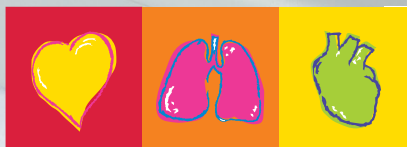

SCHILLER
The Art of Diagnostics

CARDIOVIT AT-102 G2

The new generation ECG
CARDIOVIT AT-102 G2:

- ❖ 32000 Hz sampling frequency
- ❖ Robust and reliable for hospitals and clinics
- ❖ Fast bidirectional Wi-Fi communication with HIS/EMR
- ❖ PDF report export to Pen-drive
- ❖ Best signal quality and diagnostic algorithms for adult and paediatric patients
- ❖ Spirometry and basic exercise ECG optionally available



STATE-OF-THE-ART TECHNOLOGY AND CONNECTIVITY

The CARDIOVIT AT-102 G2 features an 8" colour display with touch function keys. Moreover, its sealed alphanumeric keyboard ensures easy cleaning, making it highly suitable for the daily clinical routine. Smart Battery technology guarantees more than 8 hours of ECG recording.

Customizable printouts are available on the internal high resolution A4 thermal printer.

HIGH EFFICIENCY IN ANY WORKFLOW

ECG acquisition has never been easier and faster.



EASY 1-2-3 STEPS

- 1: Enter or download patient data
- 2: Control signal quality and lead reversal, and perform ECG
- 3: Review ECG before storing, exporting or printing

THE CARDIOVIT AT-102 G2 DIAGNOSTIC TOOLS

- ❖ Resting ECG with the SCHILLER ECG analysis program ETM as well as ETM Sport, the automatic interpretation for athletes.
- ❖ Resting rhythm recording up to 10 minutes.
- ❖ Optimal signal quality for paediatric ECGs.
- ❖ Basic exercise ECG (optional) with a large selection of protocols for treadmills and ergometers.
- ❖ SCHILLER's Culprit Coronary Artery Algorithm (CCAA, optional) locates the obstruction site in the coronary artery and thereby the size of the cardiac area at risk.

Internal high resolution
A4 thermal printer



Touch function keys

Sealed alphanumeric keyboard,
easy to clean and highly hygienic

Navigation keys allow high
efficiency in any workflow


ADD SPIROMETRY TO YOUR CARDIOVIT AT-102 G2

With the SpiroScout SP plus, the CARDIOVIT AT-102 G2 also becomes a spirometer. Thanks to GANSHORN's ultrasound technology, spirometry recordings are highly accurate and the sensor is maintenance- and calibrationfree. The disposable ScoutTubes ensure maximum hygiene and patient safety. ScoutTubes are CO₂ neutral.

The spirometry option offers FVC, SVC and MVV with a large selection of reference standards, pre and post tests, as well as interpretation of the FVC results.

The device provides feedback on the procedure's quality during the test and after each trial.

Complete spirometry results can be printed on thermal paper as well as exported to EMR/HIS systems.

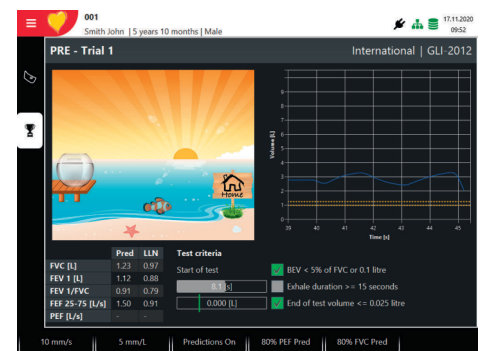

An incentive screen is included for trials with children.



SCHILLER, founded in 1974 by Alfred E. Schiller is a successful group with 31 subsidiaries and a global sales network. Today, SCHILLER is a world-leading manufacturer and supplier of devices for cardiopulmonary diagnostics, defibrillation and patient monitoring as well as software solutions for the medical industry.

Over the past 25+ years, Schiller India has established itself as a leader in the medical technology sector. With 450+ employees, a state-of-the-art ISO 13485-certified production centre in Puducherry, and multiple R&D centres across the globe, Schiller India makes advanced healthcare equipment accessible through a network of 100 sales and service dealers across more than 45 locations. Our product range includes Critical Care, Anaesthesia, Emergency Care, Cardiology, Respiratory Diagnostics, Radiology, and Robotics.



Scan here
to know more

Part No. B3.100017 01/26 Ver. 1.3

All registered trademarks acknowledged.

SCHILLER
The Art of Diagnostics

Swiss H.Q.: SCHILLER AG, Altgasse 68, P. O. Box 1052, CH -6341 Baar, Switzerland

Indian H.O.: Schiller Healthcare India Pvt. Ltd., Advance House, Makwana Road, Andheri (East), Mumbai-59. Tel.: + 91 - 22 - 61523333

No. 15/5 & 15/6, Vazhuthavur Road, Kurumbapet, Puducherry - 605 009, India

Regional Offices: Mumbai: 022-61523333 Email: sales.west@schillerindia.com Delhi: 011-41062067 / 09312432205 Email: sales.north@schillerindia.com Kolkata: 033-23593102 / 033-23593103 Email: sales.east@schillerindia.com Chennai: 044-28232648 / 044-28311021 Email: sales.south@schillerindia.com Ahmedabad: 079- 35337949 / 079- 35337950 Email: sales.west@schillerindia.com Puducherry: 0413- 220 2680 / 0413- 229 2940 Email: sales.south@schillerindia.com

CIN : U33110MH1997PTC111307

sales@schillerindia.com, help@schillerindia.com www.schillerindia.com 1800 209 8998

Follow us on : /SchillerHealthcareIndia /schiller_healthcare_india /schillerhealthcareindia /schiller_healthcare_india @schillerindia1074