


# medilogFD

The 12-lead Holter ECG for a spatial heart health analysis

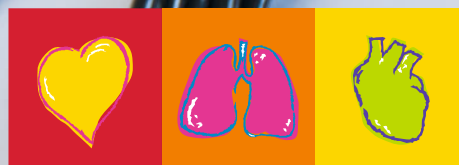

**SCHILLER**  
The Art of Diagnostics

# medilogFD

The Holter recorder that complements your medilog portfolio:

- 12-lead Holter ECG allows assessing the full spectrum of cardiac ECG conditions
- Instantaneous atrial fibrillation / atrial flutter detection based on true P-wave analysis
- Optional SpO<sub>2</sub> sensor for additional respiratory information



Easy-to-clean carrying case

## 12-LEAD HOLTER ECG FOR A SPATIAL ANALYSIS

The 12-lead Holter ECG is preferred for detailed and spatial heart analysis. Physicians, cardiologists, and electrophysiologists find great value in industry-leading technology:

- Regional assessment of myocardial ischemia. Distinguish anterior, septal, lateral, inferior region in the dynamic setting of everyday life.
- Identify the origin of ectopic beats, whether they are from the left or right ventricle.

Superior in identifying abnormal heartbeats that start in the right or left atrium.

## MADE TO SUPPRESS ARTIFACTS

The medilogFD Holter recorder features a superior resting ECG grade 128,000Hz sampling rate on 12 leads. This is used for real-time P-wave analysis, sophisticated artifact suppression, and motion detection, enabling a speedy evaluation.

## HOLISTIC APPROACH

Instantaneous detection of atrial fibrillation onset, analysis of respiratory episodes, and high resolution HRV are just some of the exceptional features of the medilogFD.

## DUAL-BATTERY CONCEPT

The high-quality internal battery is recharged via USB or wirelessly on an optional charging station. Moreover, easily available standard AAA batteries can be used. This allows patients to be screened for up to five days on a single charge.

## COMPLETES THE MEDILOG PORTFOLIO

medilogFD completes your medilog line-up with medilogAR (3-lead Holter ECG, up to 14 days runtime) and medilog DARWIN2 (high-end ECG analysis software).

High-contrast colour display for easy operation. Live ECG preview: 3 channels shown at a time, all 12 can be scrolled through

Multi-factor signal check based on:

- Impedance
- Amplitude
- Amount of noise

Microphone to record patient ID. In an anonymous start, a voice recording is a must to prevent patient mix-up

Shock-proof thanks to its robust design

Large and easy-to-operate buttons

Visible and tactile position markers for easy patient cable connection

Eco-friendly: runs on an internal rechargeable battery, or, if needed, on normal AAA batteries for fast turnaround times

Data exchange via USB with robust state-of-the-art USB-C connector

Removable patient cable with 10 wires, tested for more than 4000 connection cycles



Wireless SpO<sub>2</sub> sensor, attached for the night

## SCREEN FOR RESPIRATION EPISODES

Thanks to ECG-derived respiration recording (EDR), the medilogFD can screen for potential respiratory episodes during sleep. The optional SpO<sub>2</sub> sensor, connected via Bluetooth, allows for additional respiratory information.



## About SCHILLER

SCHILLER A. G., a company based in Switzerland, was founded by Alfred Schiller in 1974. Over the last 50 years, it has evolved into a successful Global Group with 31 subsidiaries and dealers across 30 countries. Today, SCHILLER is a world-leading manufacturer and supplier of devices for cardiopulmonary diagnostics, defibrillation and patient monitoring as well as software solutions for the medical industry.

Schiller in India was founded as a joint venture in 1997, and since then has established itself as a leader in the medical technology sector. With 450+ employees, a state-of-the-art ISO 13485-certified production centre in Puducherry, multiple R&D centres, and 17 offices across India; Schiller makes advanced healthcare equipment accessible through a network of 100 sales and service dealers across more than 45 locations. Our product range includes Critical Care, Anaesthesia, Emergency Care, Cardiology, Respiratory Diagnostics, Radiology, and Robotics.

Scan here  
to know more



Part No. B3.100032 07/2025 Ver. 1

All registered trademarks acknowledged.

**SCHILLER**  
The Art of Diagnostics

Swiss H.Q.: SCHILLER AG, Altgasse 68, P. O. Box 1052, CH -6341 Baar, Switzerland

Indian H.O.: Schiller Healthcare India Pvt. Ltd., Advance House, Makwana Road, Andheri (East), Mumbai-59. Tel: + 91 - 22 - 61523333

🏠 : No. 15/5 & 15/6, Vazhuthavur Road, Kurumbapet, Puducherry - 605 009, India

Regional Offices: Mumbai: 022-61523333 Email: sales.west@schillerindia.com Delhi: 011-41062067 / 09312432205 Email: sales.north@schillerindia.com Kolkata: 033-23593102 / 033-23593103 Email: sales.east@schillerindia.com Chennai: 044-28232648 / 044-28311021 Email: sales.south@schillerindia.com Ahmedabad: 079- 35337949 / 079- 35337950 Email: sales.west2@schillerindia.com Puducherry: 0413- 2202680 / 0413- 229 2940 Email: sales.south@schillerindia.com

CIN : U33110MH1997PTC111307

✉ sales@schillerindia.com, help@schillerindia.com 🌐 www.schillerindia.com 📞 1800 209 8998

Follow us on : [f](#) /SchillerHealthcareIndia [i](#) /schiller\_healthcare\_india [in](#) /schillerhealthcareindia [c](#) /schiller\_healthcare\_india [v](#) @schillerindia1074