


Comprehensive monitoring
Volumetric capnography
High-flow oxygen therapy
Spontaneous Breathing Test

High-end mechanical ventilator that
provides invasive and non-invasive
ventilation, adaptable to any
patient category

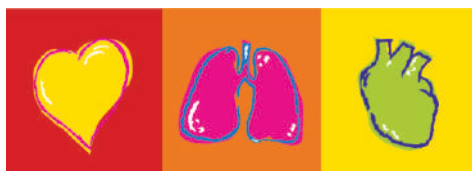

SCHILLER
The Art of Diagnostics

GraphNet advance

GraphNet advance is a high-end mechanical ventilator providing both invasive and non-invasive ventilation, tailored to suit any patient category. It delivers precise and reliable respiratory support for a wide range of clinical applications.



COMPACT DESIGN

- ❖ Equipped with a built-in 12.1" LED touch screen
- ❖ Offers a clear and comprehensive display of vital sign variables
- ❖ Ensures easy monitoring of patient data at a glance
- ❖ User-friendly and intuitive programming with high visibility alarm indicator
- ❖ Has a fully integrated front panel without gaps and is easy to clean



INTEGRATED MONITORING

- ❖ Comprehensive patient monitoring is built into the system, offering real-time feedback on critical parameters like lung compliance, resistance, and tidal volume, enabling clinicians to make informed decisions promptly.

ADAPTIVE VENTILATORY ASSISTANCE (AVA)

- ❖ This advanced feature allows for automatic adjustment of the ventilatory parameters in response to the patient's needs, promoting lung-protective strategies and reducing the risk of complications.



- ❖ Supports ventilation for adult, paediatric, and neonatal patients
- ❖ Offers high-flow oxygen therapy
- ❖ Features automatic leak compensation
- ❖ Includes integrated volumetric capnography software
- ❖ Suitable for intra-hospital transportation
- ❖ Equipped with Spontaneous Breathing Test (SBT)
- ❖ Low maintenance cost
- ❖ Built-in battery with up to 2.5 hours of capacity



Volumetric Capnography

Integrated monitoring of exhaled CO_2 and display of volumetric curves over time, as well as measurements of physiological dead space.



Spontaneous Breathing Test

This new tool allows a spontaneous breathing test to be performed before weaning, monitored in real time and completely safe with dedicated alarms.



High-flow Oxygen Therapy

The addition of this method expands the use of the ventilator to the stages that are prior to and after the mechanical ventilation itself, by providing non-invasive assistance to oxygenation.



AVA (Adaptive Ventilatory Assistance)

The AVA mode is an adaptive & closed-loop with optimal control that guarantees minute ventilation of the patient, generating a ventilatory pattern with minimal respiratory effort.





About SCHILLER

SCHILLER A. G., a company based in Switzerland, was founded by Alfred Schiller in 1974. Over the last 50 years, it has evolved into a successful Global Group with 31 subsidiaries and dealers across 30 countries. Today, SCHILLER is a world-leading manufacturer and supplier of devices for cardiopulmonary diagnostics, defibrillation and patient monitoring as well as software solutions for the medical industry.

Schiller in India was founded as a joint venture in 1997, and since then has established itself as a leader in the medical technology sector. With 450+ employees, a state-of-the-art ISO 13485-certified production centre in Puducherry, multiple R&D centres, and 17 offices across India; Schiller makes advanced healthcare equipment accessible through a network of 100 sales and service dealers across more than 45 locations. Our product range includes Critical Care, Anaesthesia, Emergency Care, Cardiology, Respiratory Diagnostics, Radiology, and Robotics.

Scan here
to know more:



Part no.: B11.100013  07/26 Ver. no.: 1.2

All registered trademarks acknowledged. Accessories shown are not part of standard configuration

SCHILLER
The Art of Diagnostics

Swiss H.Q.: SCHILLER AG, Altgasse 68, P. O. Box 1052, CH -6341 Baar, Switzerland

Indian H.O.: Schiller Healthcare India Pvt. Ltd., Advance House, Makwana Road, Andheri (East), Mumbai-59. Tel : + 91 - 22 - 61523333

 : No. 15/5 & 15/6, Vazhuthavur Road, Kurumbapet, Puducherry - 605 009, India

Regional Offices: Mumbai: 022-61523333 Email: sales.west@schillerindia.com Delhi: 011-41062067 / 09312432205 Email: sales.north@schillerindia.com Kolkata: 033-23593102 / 033-23593103 Email: sales.east@schillerindia.com Chennai: 044-28232648 / 044-28311021 Email: sales.south@schillerindia.com Ahmedabad: 079- 35337949 / 079- 35337950 Email: sales.west2@schillerindia.com Puducherry: 0413- 220 2680 / 0413- 229 2940 Email: sales.south@schillerindia.com

CIN : U33110MH1997PTC111307

 sales@schillerindia.com, help@schillerindia.com  www.schillerindia.com  1800 209 8998

Follow us on :  /SchillerHealthcareIndia  /schiller_healthcare_india  /schillerhealthcareindia  /schiller_healthcare_india  @schillerindia1074